




TooToo MediTech Co., Ltd.



Wi-Fi
(optional)

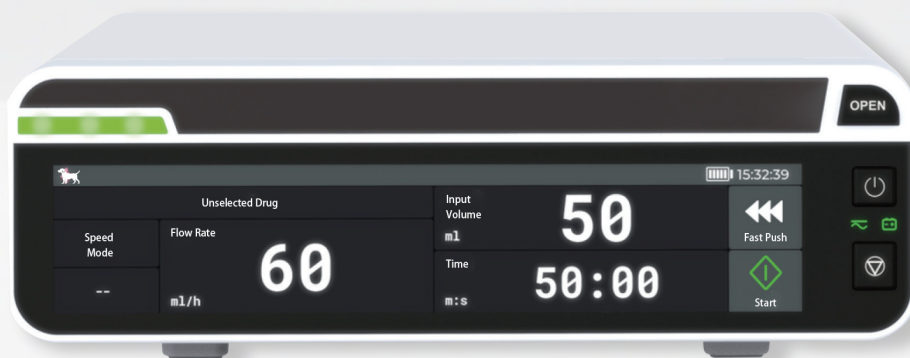

VP50

Veterinary Infusion System

Smart touch, Safe medication



7" full-size touchscreen
bringing a brand new user experience



See-through window design
easy to clean and observe syringe status



Intuitive UI design
80% operation can be done
on the main interface



Premium material
supports up to 49 disinfectants

Address: 3A, Phase II, Smart Home, Baolong Street,
Longgang District, Shenzhen, Guangdong, China
Tel: +86-755-28306385 Email: info@newtootoo.com
Web: www.newtootoo.com



Touch the flexibility

Extraordinary touch experience


Infusion System Series features a full-size color touchscreen that comes with an intuitive UI to provide an extraordinary interactive experience as never before.

- 7-inch HD capacitive touchscreen makes the interaction as smooth as using smart devices in your daily life.
- The intuitive UI facilitates an easy operation by presenting critical parameters with an elegant layout.



All in one


Infusion System sSP and sIP are versatile enough to integrate all functions to satisfy a variety of clinical purposes.

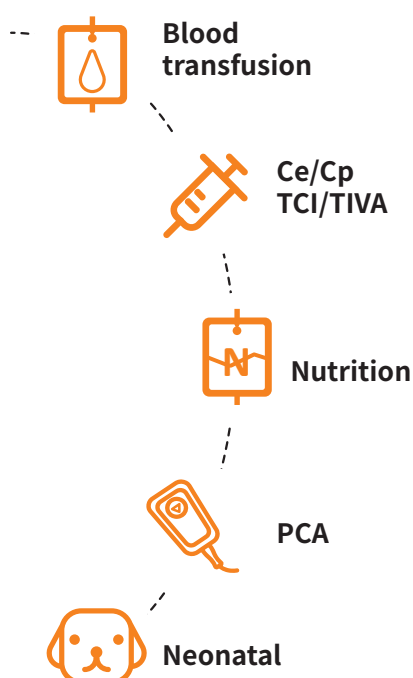




Dynamic Pressure System(DPS)

Speedometer style indicator with numerical pressure to monitor the in-line pressure trend just with a glance.

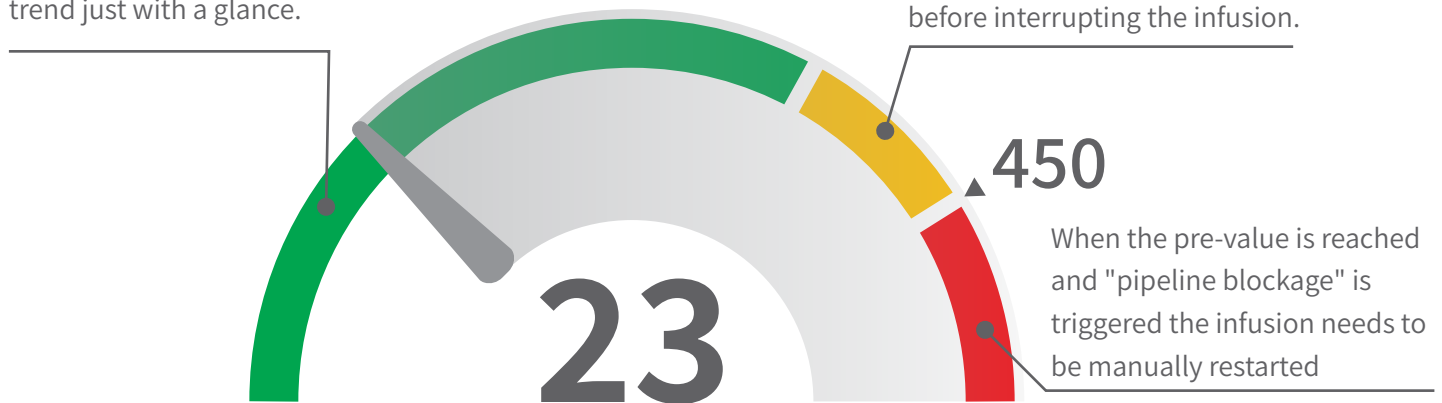


Guided Alarm Management


The highly visible alarm management enables users to fast locate and solve the problem.

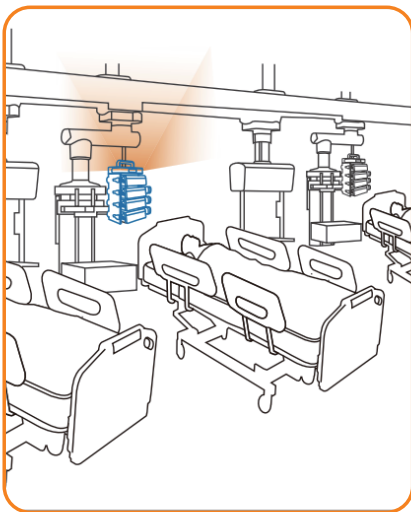

Fast location of the bedside alarm with 270° visible alarm light on the Infusion System from a distance.

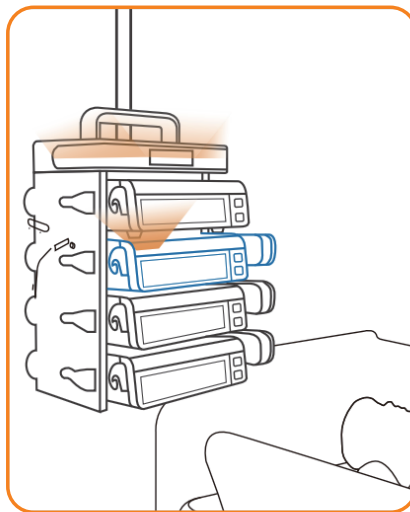

Fast location of the pump with an audible and visual alarm on Infusion System at bedside.

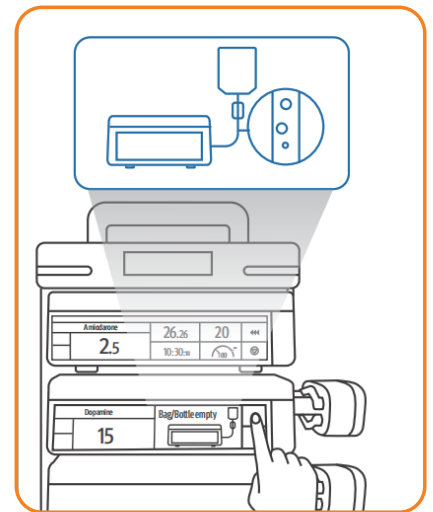

Fast location of the problem with an instructional animation displayed on the pump to aid problem resolution



Managing drugs safely throughout the infusion

Color-coding system

The color coding of drug name assists users to easily select and verify the correct drug.

. SafeDosen Info software enables programming infusion parameters automatically to enhance efficiency.



SafeDose™ DERS

The SafeDose™ DERS helps to reduce drug errors by preventing users from exceeding pre-set hard limits, and by requiring confirmation to override soft limits.



Program the customized drug information with color and limits.



Select the drugs easily with clear categorization and search function.



A clear alert message is displayed if a limit is reached.



Dedicated to safety

Innovative technologies ensures safe medication


SmartRapid

SmartRapid ensures timely infusion by significantly shortening the start-up time, from turning on the pump to delivering the first drop of medication.



Precise and consistent performance

- Rivalled accuracy: $sSP \leq +1.8\%$; $sIP \leq +5\%$
- Flow rate initiates as precise as 0.01ml/h
- Compatible syringes as small as 1ml

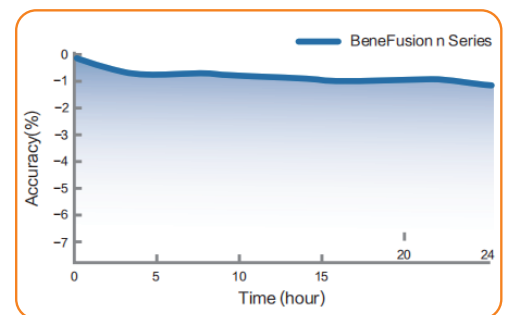

Integrated central monitoring

Infusion System offers one-stop monitoring of all pets' vital sign and infusion treatment details, providing comprehensive information for clinical workers to improve the quality of pet care.



Closed-loop medication management

Assists clinical staff to reduce medication errors with barcode scanner by aligning prescription automatically, enhancing the work-flow efficiency.

